


Control, Durability, and Versatility – All in One.

Traxcess® Guidewire, the first neuro-guidewire with hybrid nitinol/stainless steel technology, is uniquely designed to reach challenging distal anatomies. It combines exceptional tip softness and flexibility with precise torque response, to help you get to your destination.




Control

- Hybrid nitinol/stainless steel construction offers a soft, atraumatic tip while maximizing push/pull control and torque response*
- Lubricious hydrophilic coating allows smooth advancement of the guidewire within the vasculature and microcatheter*
- Excellent tip shapability and shape retention design promotes controlled steerability*

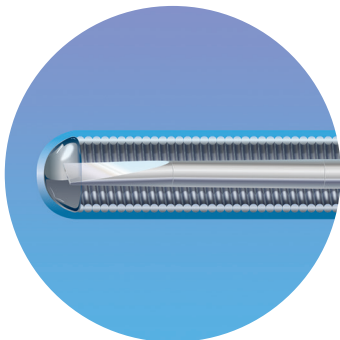


Durability

- Proprietary hydrophilic coating offers exceptional integrity and durability*
- Kink-resistant distal nitinol core maintains its shape throughout the procedure, even through the most tortuous vessels*

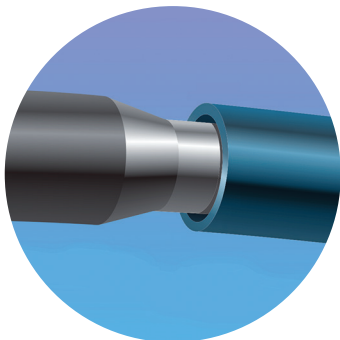

Versatility

- Innovative Traxcess® Docking Wire simplifies catheter exchange, intra-procedurally, reducing time and cost

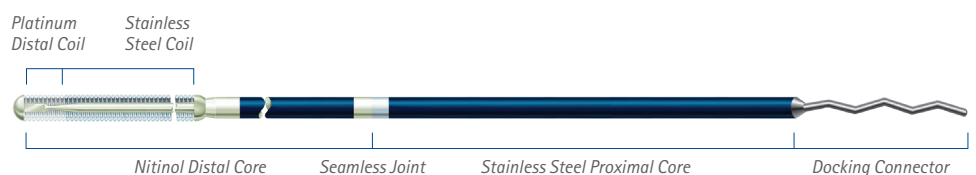

Traxcess[®] 14 Guidewire

Product Name	Description	Product Code	Total Length	Hydrophilic Coating Length	Radiopaque Tip Length	Distal OD	Proximal OD
Traxcess [®] 14	Soft	GW1420040	200 cm	40 cm	3 cm	0.012 in (0.30 mm)	0.014 in (0.36 mm)
Traxcess [®] 14 SELECT	Soft	GW1420040S	200 cm	97 cm	3 cm	0.012 in (0.30 mm)	0.014 in (0.36 mm)
Traxcess [®] 14 EX	Support	GW1420040X	200 cm	40 cm	6 cm	0.012 in (0.30 mm)	0.014 in (0.36 mm)

Straight, shapeable tip. Shapeable tip length: 1.4 cm.
Packed 1 guidewire per box. Includes shaping mandrel, insertion tool, and torque device.

Added features of the new Traxcess[®] 14 SELECT Guidewire

- Distal 3 cm tip redesigned for improved support and better vessel selectivity
- Longer hydrophilic coating for enhanced overall navigation and reduced friction with compatible catheters
- Optimized accessories: longer wire shaping mandrel, torque device with improved grip, and longer introducer

Traxcess[®] Docking Wire

Product Name	Product Code	Total Length	Outer Diameter
Traxcess [®] Docking Wire	GW14100EX	115 cm	0.014 in (0.36 mm)

Packed 1 docking wire per box.

- Traxcess[®] Docking Wire extends the length of Traxcess[®] 14 Guidewires into exchange length (313 cm) without the need for removing the guidewire, which is then easily returned to the original 200 cm by removing the docking wire

For more information or a product demonstration, contact your local MicroVention representative.



INDICATIONS FOR USE:

The Traxcess[®] 14 Guidewire is intended for general intravascular use, including the neuro and peripheral vasculature. The guidewire can be steered to facilitate the selective placement of diagnostic or therapeutic catheters. The device is not intended for use in coronary arteries. There are no known contraindications. Refer to Instructions for Use for additional information.

MicroVention, Inc.

Worldwide Headquarters

1311 Valencia Avenue
Tustin, CA 92780 USA

MicroVention UK Limited

MicroVention Europe, S.A.R.L.

MicroVention Deutschland GmbH

Web

PH +1.714.247.8000

PH +44 (0) 191 258 6777

PH +33 (1) 39 21 77 46

PH +49 211 210 798-0

microvention.com